



Millimeter-wave therapy instrument Using guide

Contents

Contents	- 3 -
I. Millimeter wave therapy	- 5 -
1. Health maintenance, and health recuperation:	- 5 -
2. Diabetic main site:.....	- 6 -
3. Diabetic foot:	- 7 -
4. Diabetic skin ulcer:	- 7 -
5. Diabetic neuropathy, Vascular lesions :	- 8 -
6. Diabetic nephropathy:.....	- 8 -
7. Diabetic retinopathy:	- 9 -
8. Diabetic thirst:	- 10 -
9. Diabetic hunger:	- 11 -
10. Tired and weak, few words, diarrhea:.....	- 12 -
11. Lower extremity weakness, cramps, cold and weakness:.....	- 13 -
12. Prostatitis, Hyperplasia of prostate gland:	- 14 -
13. Mammitis, Hyperplasia of mammary glands:.....	- 14 -
14. Trauma and Wound healing:	- 14 -
15. Acute soft tissue injury:	- 15 -
16. Suppurative inflammation:	- 15 -
17. Erysipelas:	- 15 -
18. Acute Angialymphitis, Adenolymphitis:	- 15 -
19. Cervical spondylosis:	- 15 -
20. Scapulohumeral periartthritis:	- 16 -
21. Tennis elbow, Tendonitis:	- 16 -
22. Lumbar sprain, Lumbar muscle strain:.....	- 16 -
23. lumbar disc protrusion:	- 16 -
24. Osteoarthritis:	- 16 -
25. ONFH (Osteonecrosis of the Femeral Head) :	- 16 -
26. Temporomandibular joint(Pic 35) disorder:.....	- 16 -
27. Allergic rhinitis, Nasal polyps, Inferior turbinate hypertrophy:.....	- 17 -
28. Toothache:	- 17 -
29. Sore throat:.....	- 17 -
30. Headaches, trigeminal neuralgia, facial palsy, hemifacial spasm:	- 17 -
31. Bronchitis, bronchial asthma:	- 18 -
32. Pain in the stomach, gastritis, peptic ulcer:.....	- 18 -
33. Hypertension:.....	- 19 -


34. Cerebral infarction:	- 19 -
35. Coronary heart disease, angina pectoris:.....	- 20 -
36. Tonsillitis:	- 20 -
37. Nephropathy:	- 21 -
38. Urinary tract infection:	- 22 -
39. Constipation:.....	- 22 -
40. Insomnia, neurasthenia.....	- 23 -
41. Dysmenorrhea:	- 23 -
42. Pelvic inflammation、Pelvic inflammatory mass、Accessory phlogistic、Senile vaginitis:	- 24 -
43. Acne:	- 24 -
44. Eczema, Herpes zoster, Herpes zoster neuralgia, Allergic dermatitis	- 24 -
II. Millimeter-wave cancer treatment.....	- 25 -
1. Cancer of head and neck:	- 25 -
2. Breast cancer:	- 26 -
3. Lung cancer:	- 26 -
4. Rectal cancer:	- 27 -
5. Bladder cancer:.....	- 28 -

I. Millimeter wave therapy

(Millimeter wave probe is big, you can cover the accupoint with it.)


1. Health maintenance, and health recuperation:

Irradiate **Zusanli(Pic 1)**, **Hegu(Pic 2)** and **Neiguan(Pic 3)** accupoint.
One or two times a day, 15 minutes each time.

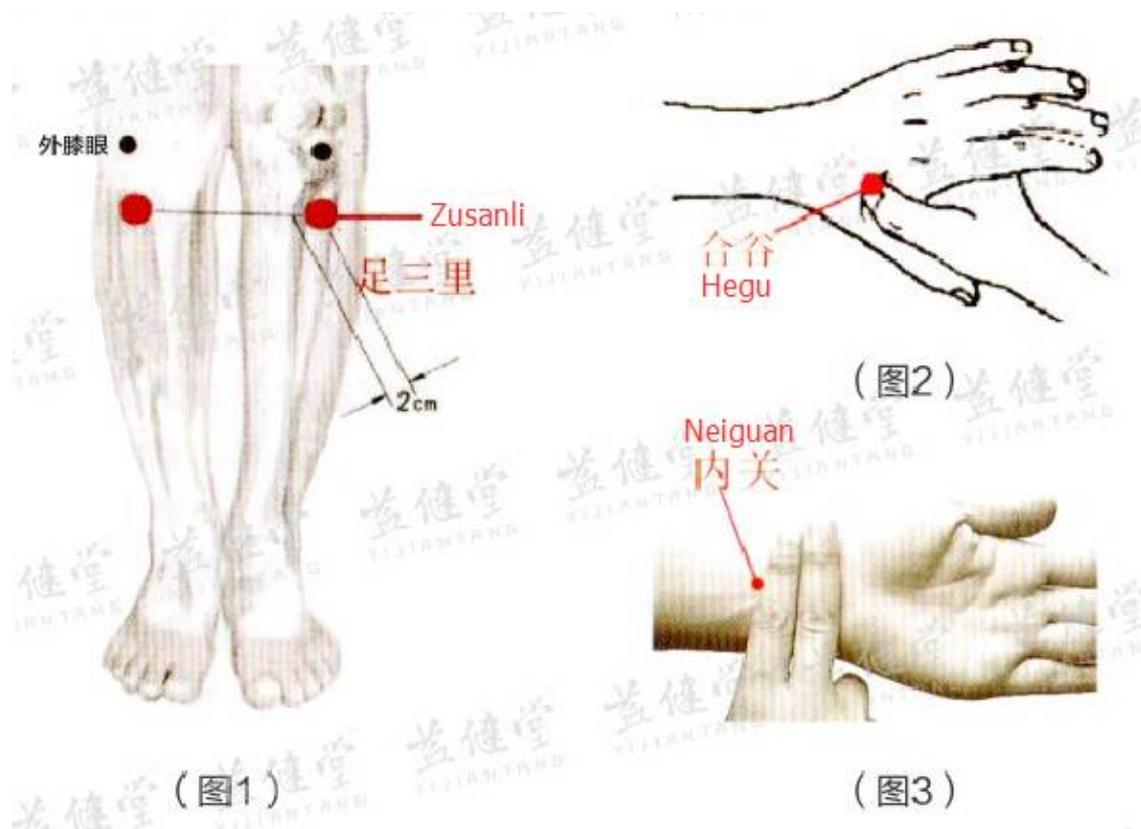

2. Diabetic main site:

Irradiate of pancreatic body surface projection area, **Diji(Pic 4)**, **Shenshu(Pic 5)**, **Sanyinjiao(Pic 6)** accupoint. One or two times a day, 30 minutes each time.

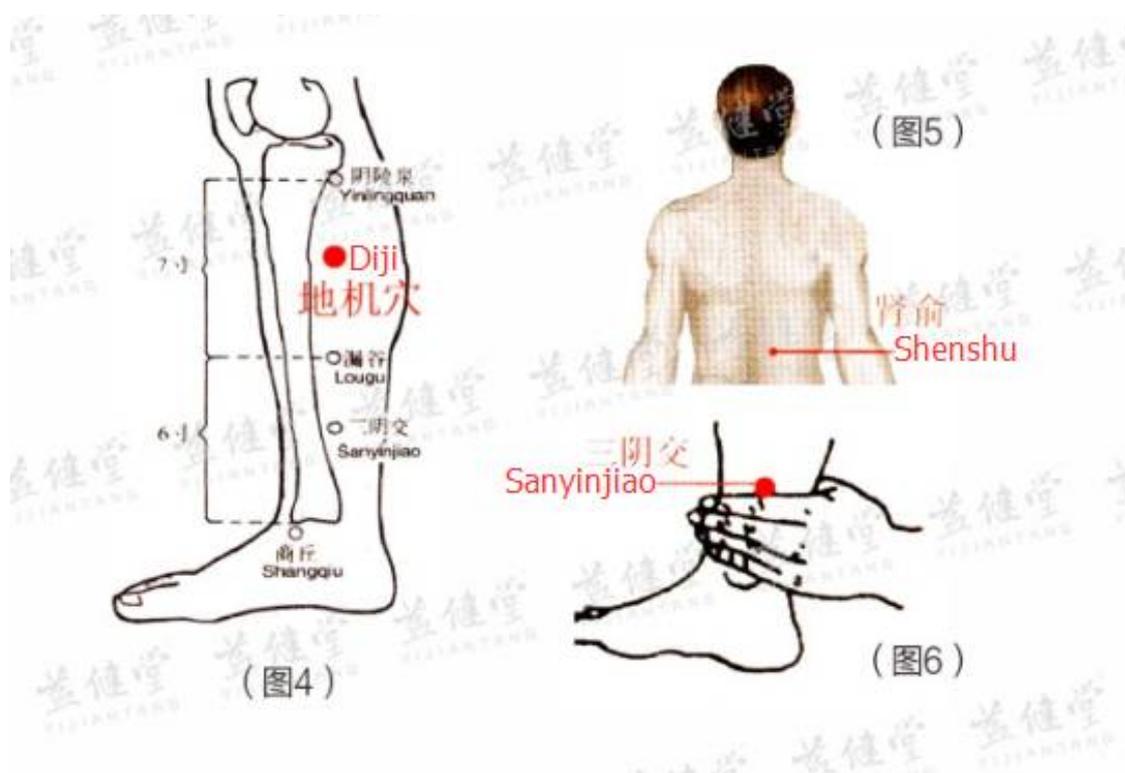


3. Diabetic foot:

Irradiate lesions area, auxiliary the **surface projection of the pancreas(Pic 7)**, **Diji(Pic 8) accupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.

4. Diabetic skin ulcer:

Irradiate lesions area, auxiliary the **surface projection of the pancreas(Pic 7)**, **Diji(Pic 8) accupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.

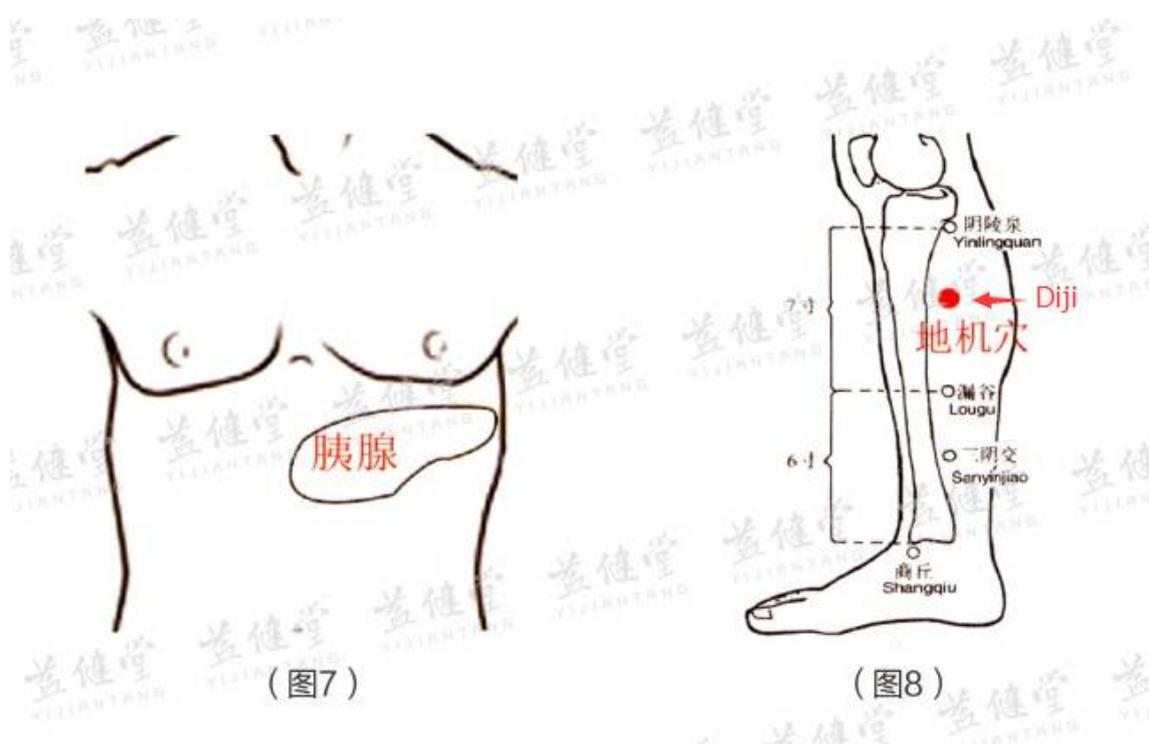


5. Diabetic neuropathy, Vascular lesions :

Irradiate pain area, auxiliary the **surface projection of the pancreas(Pic 7), Diji(Pic 8), Neiting(Pic 9) accupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.


6. Diabetic nephropathy:

Irradiate **Shenshu(Pic 10)**, auxiliary the **surface projection of the pancreas(Pic 7), Diji(Pic 8) accupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.

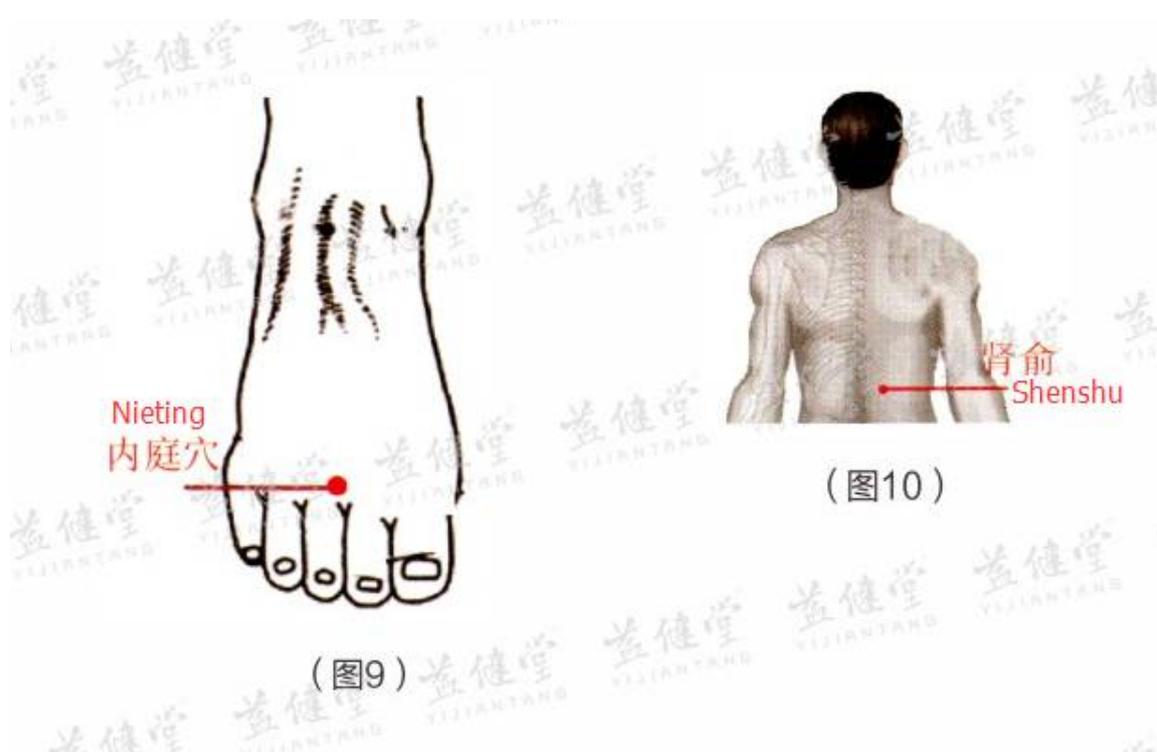


7. Diabetic retinopathy:

Irradiate the **surface projection of the pancreas(Pic 11), Diji(Pic 12), Ganshu(Pic 13)**, auxiliary **Guangming(Pic 14) accupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.

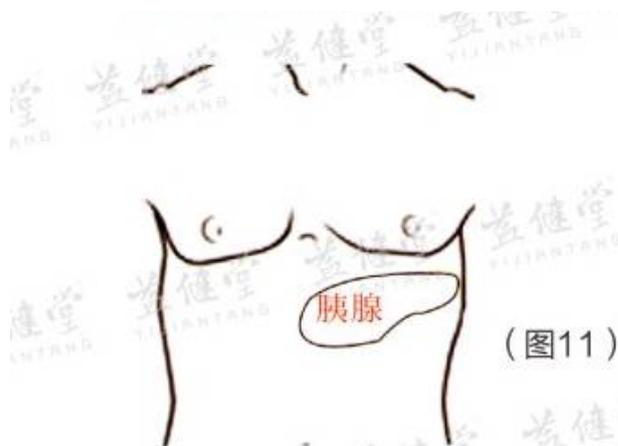

8. Diabetic thirst:

Irradiate **Yishe(Pic 15)**, **Feishu(Pic 16)**, **Chengjiang(Pic 17)**, auxiliary the **surface projection of the pancreas(Pic 11)**, **Diji(Pic 12)** accupoint. One or two times a day, 30 minutes each part. Continue treatment until recovery.

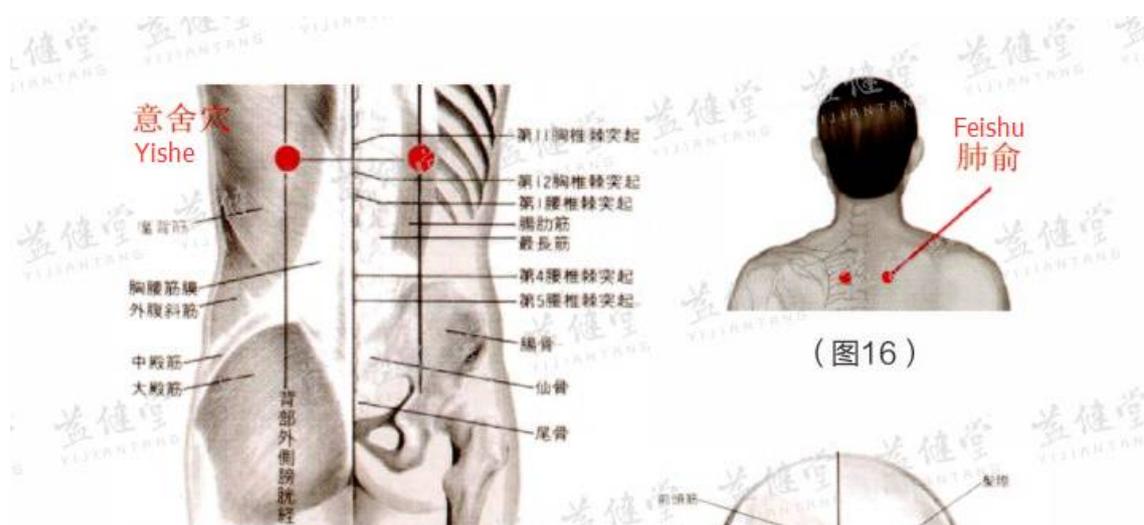


(图16)



(图17)


9. Diabetic hunger:

Irradiate **Fenglong(Pic 18), Diji(Pic 19) accupoint**, auxiliary the **surface projection of the pancreas(Pic 20)**, one or two times a day, 30 minutes each part. Continue treatment until recovery.

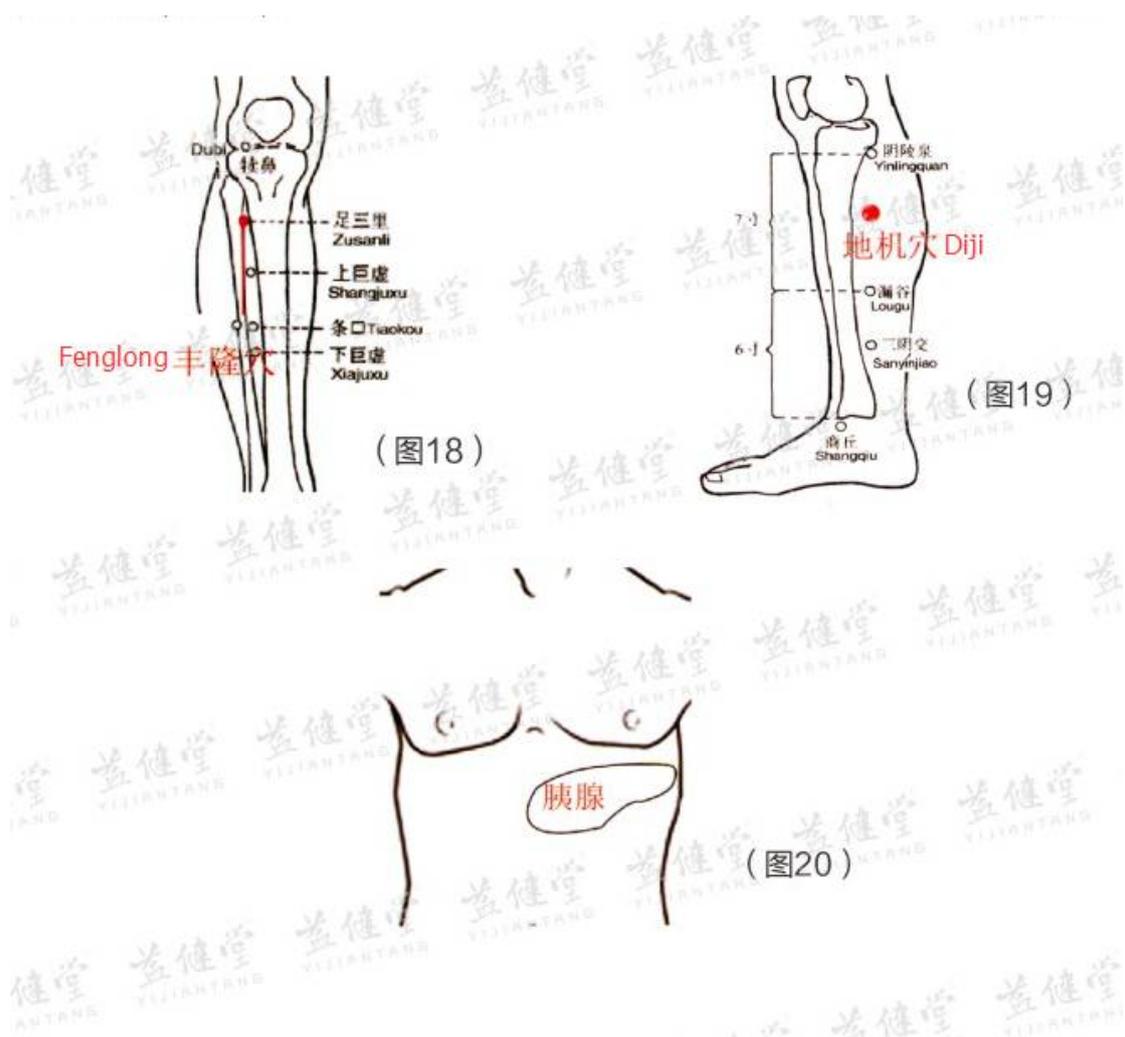


10. Tired and weak, few words, diarrhea:

Irradiate **Diji(pi19) acupoint**, auxiliary the **surface projection of the pancreas(Pic 20)**, **Weishu(Pic 21)**, **Yinlingquan(Pic 22)**, **Sanyinjiao(Pic 23) acupoint**. One or two times a day, 30 minutes each part. Continue treatment until recovery.

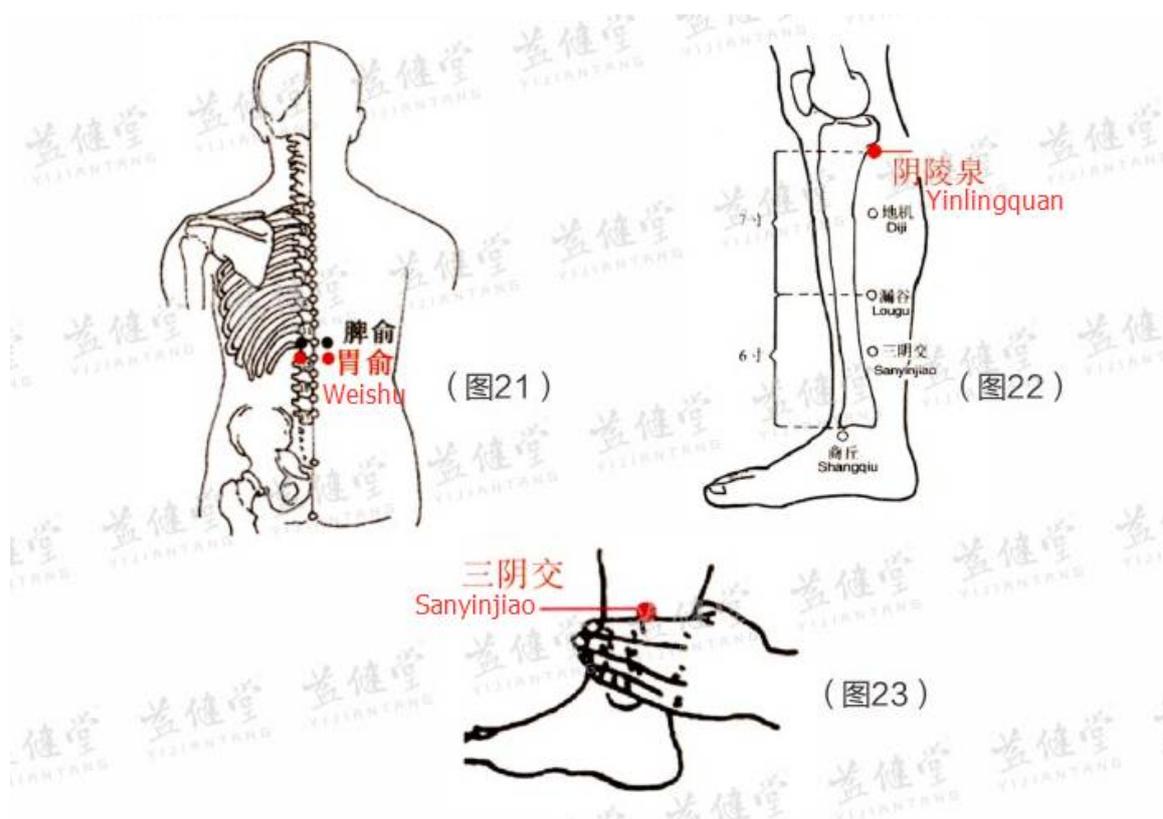

11. Lower extremity weakness, cramps, cold and weakness:

Irradiate **Xuehai(Pic 27)**, **Zusanli(Pic 28)**, **Sanyinjiao(Pic 29)** acupoint, and the affected area. One or two times a day, 30 minutes each part. Continue treatment until recovery.

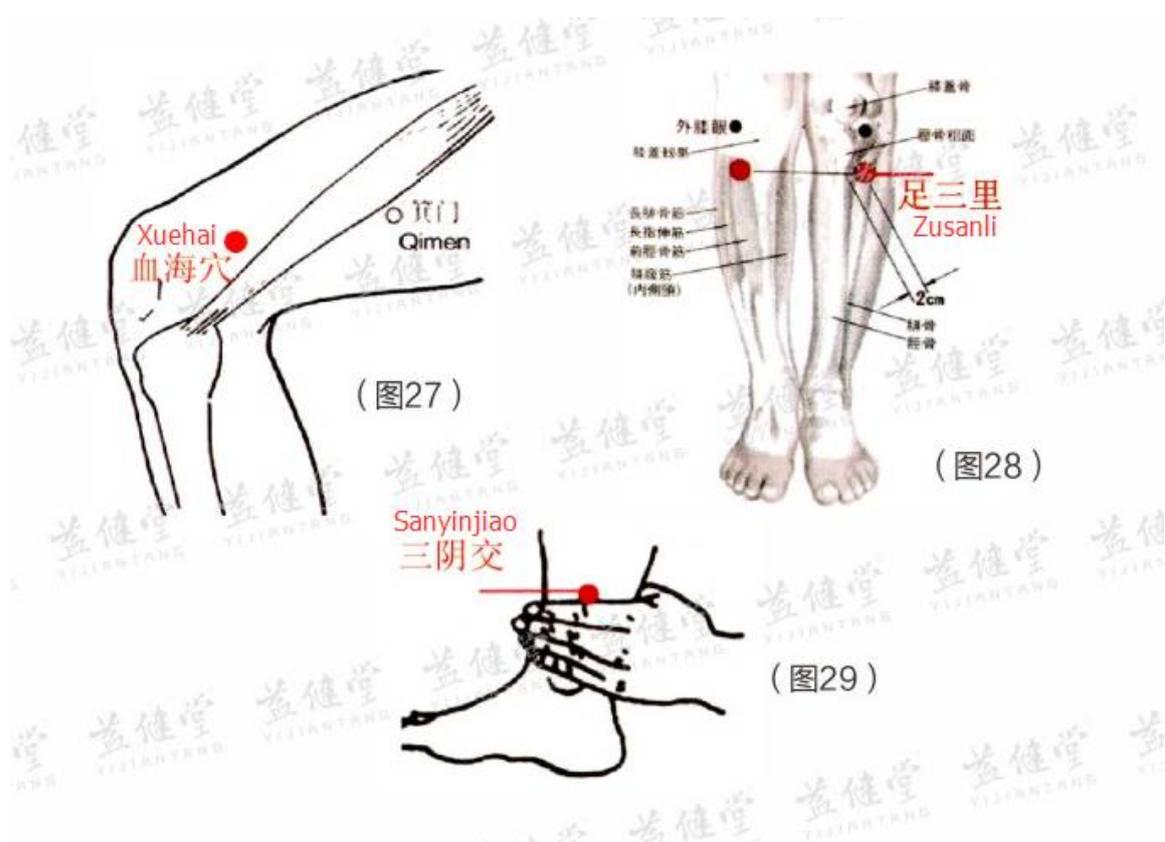

12. Prostatitis, Hyperplasia of prostate gland:


Irradiate **Huiyin(Pic 30)**. 1 or 2 times a day, 30 minutes to 60 minutes each time. Continue treatment until recovery.

13. Mammitis, Hyperplasia of mammary glands:

Irradiate lesions area, auxiliary **Rugen(Pic 31) or Danzhong(Pic 32) accupoint**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

14. Trauma and Wound healing:

Irradiate wounds and lesions area. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

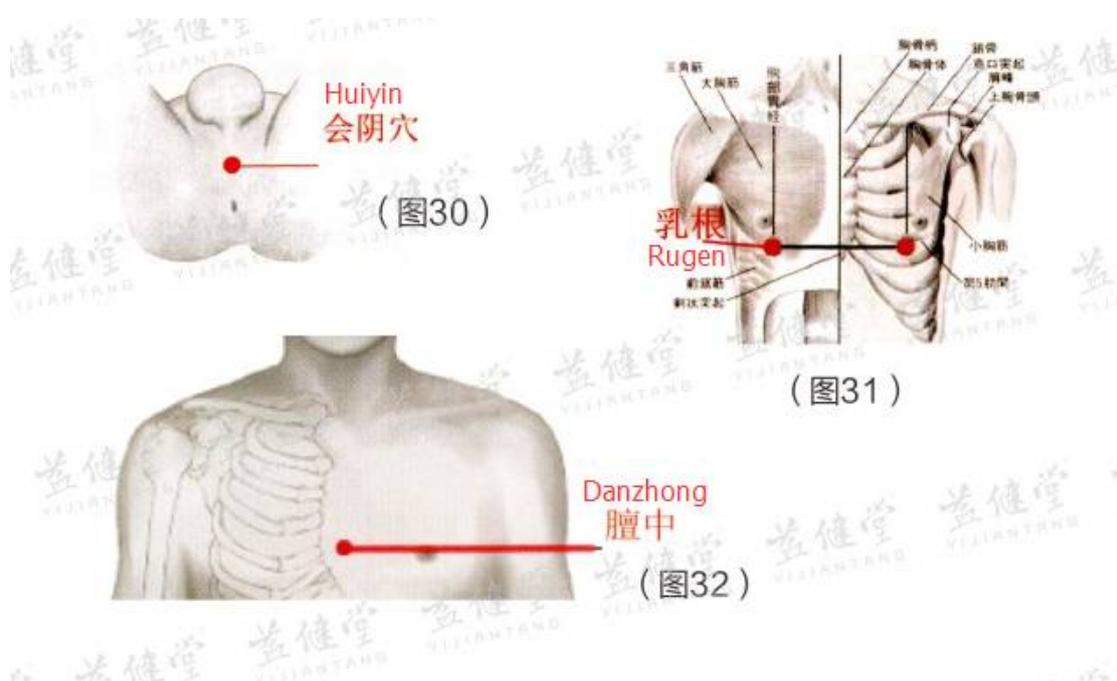

15. Acute soft tissue injury:

Irradiate injury, 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

16. Suppurative inflammation:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

17. Erysipelas:


Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

18. Acute Angialymphitis, Adenolymphitis:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

19. Cervical spondylosis:

Irradiate **Dazhui(Pic 33),Fengchi(Pic 34) accupoint** and lesions area. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

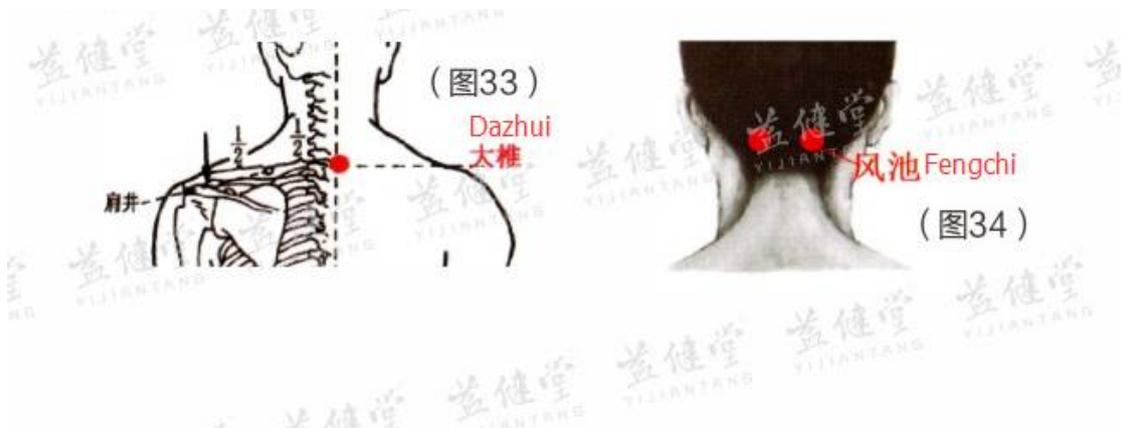

20. Scapulohumeral periarthritiis:

Irradiate pain area, 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

21. Tennis elbow, Tendonitis:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

22. Lumbar sprain, Lumbar muscle strain:

Irradiate pain area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

23. Lumbar disc protrusion:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

24. Osteoarthritis:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

25. ONFH (Osteonecrosis of the Femeral Head) :

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.

26. Temporomandibular joint(Pic 35) disorder:

Irradiate pain area. 1 or 2 times a day, 30 minutes each time.
Continue treatment until recovery.



(图35)

27. Allergic rhinitis, Nasal polyps, Inferior turbinate hypertrophy:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time.

Continue treatment until recovery.

28. Toothache:

Irradiate pain area and **Hegu(Pic 36) accupoint**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

29. Sore throat:


Irradiate lesions area, auxiliary **Tiantu(Pic 37) accupoint**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

30. Headaches, trigeminal neuralgia, facial palsy, hemifacial spasm:


Irradiate lesions area, auxiliary **Taiyang(Pic 37), Zusanli(Pic 38) accupoint**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.



(图36)

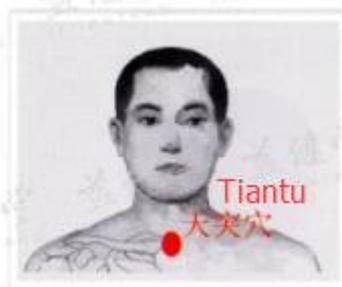

(图37)

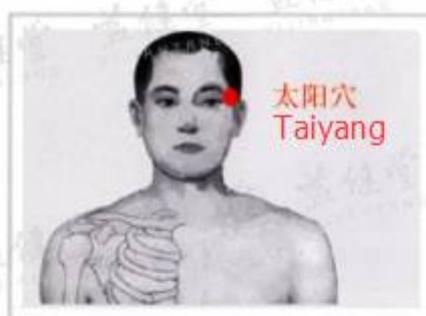

(图38)

31. Bronchitis, bronchial asthma:

Irradiate lesions area, **Dingchuan(Pic 10), Tiantu(Pic 41)** accupoint. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

32. Pain in the stomach, gastritis, peptic ulcer:


Irradiate stomach lesions area, **Zusanli(Pic 41)** accupoint. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.



(图39)



(图40)

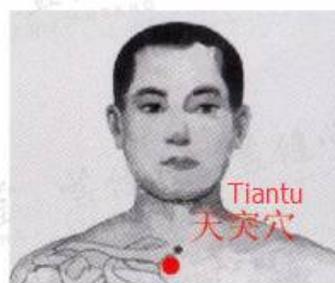


(图41)

33. Hypertension:

Irradiate carotid artery area, **Quchi(Pic 42)**, **Neiguan(43)** accupoint. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

34. Cerebral infarction:

Irradiate carotid artery area, **Quchi(Pic 42)**, **Neiguan(43)** accupoint and **right shoulder joint area(Pic 44)**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.


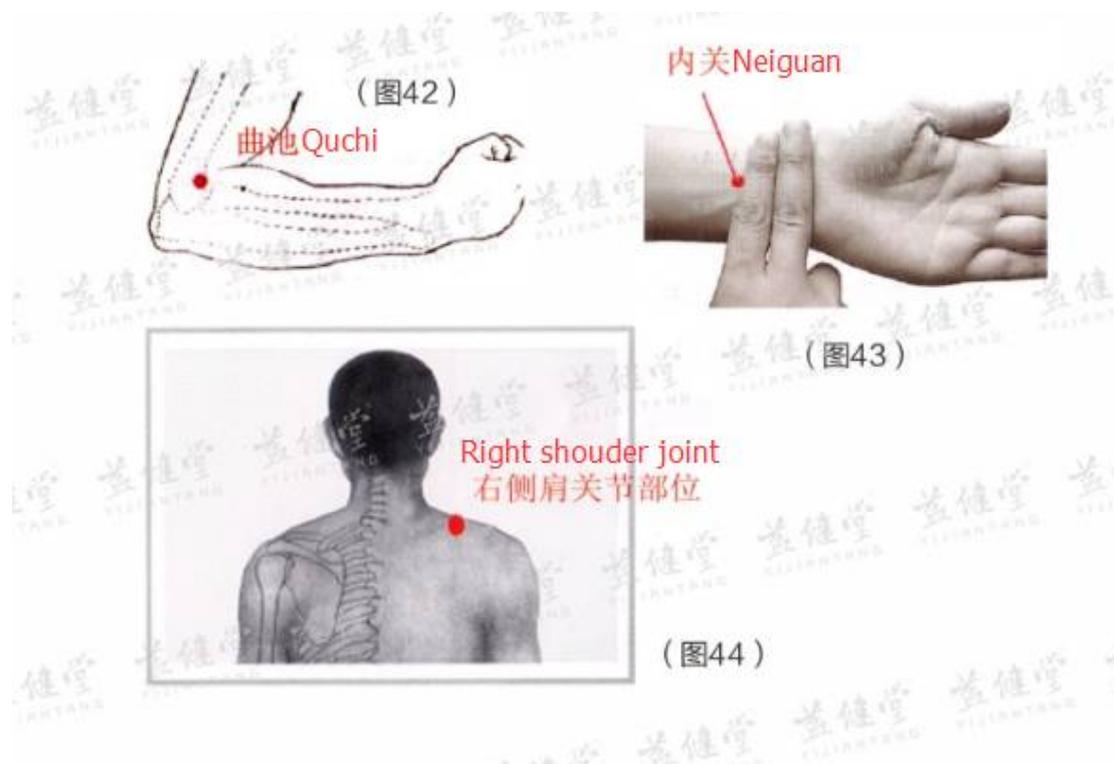


35. Coronary heart disease, angina pectoris:

Irradiate Chest region, **Neiguan(Pic 45)**, **Xinshu(46)** accupoint. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.


36. Tonsillitis:

Irradiate lesions area. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

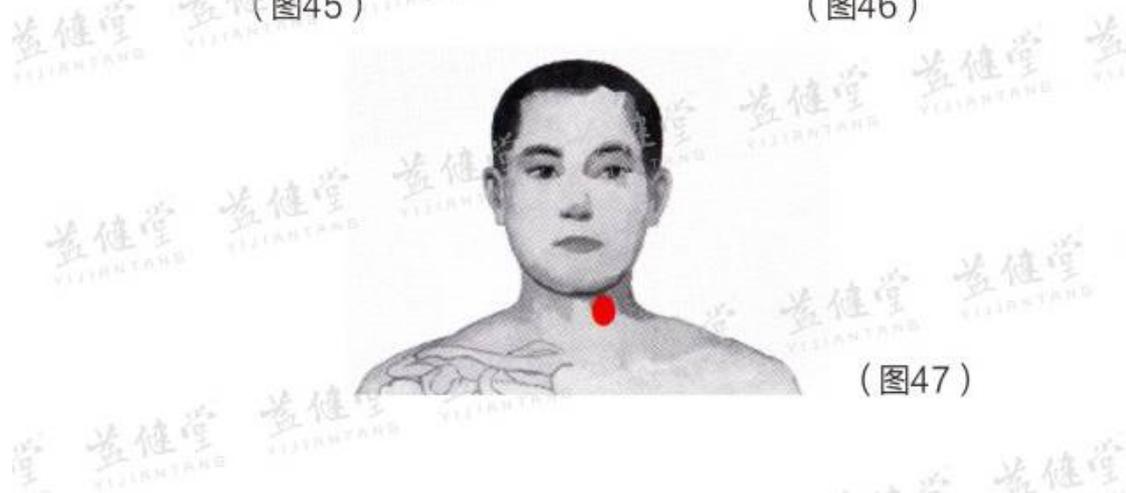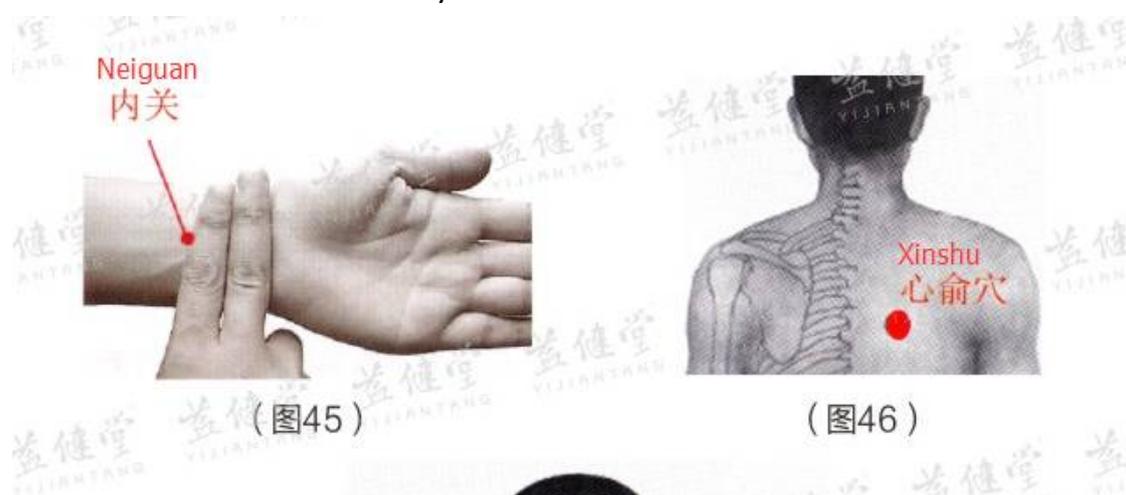

37. Nephropathy:

Irradiate **Shenshu**(Pic 45), **Zusanli**(46) accupoint. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

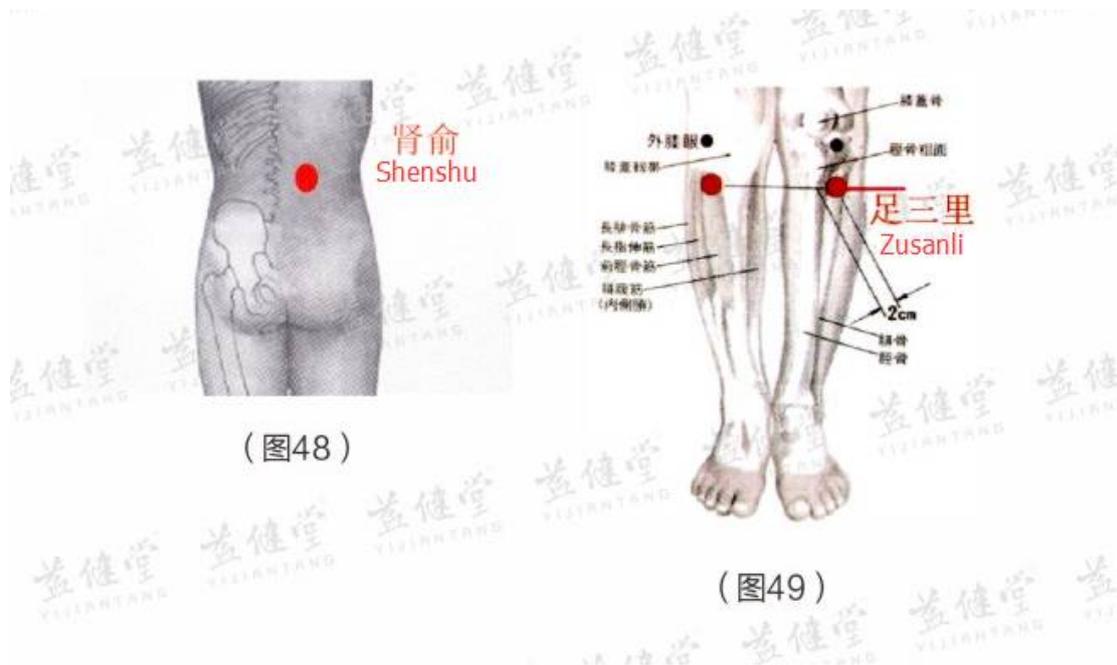


38. Urinary tract infection:

Irradiate both **underbelly, Yinlingquan acupoint (Pic 50)** , 1 or 2 times a day, 30-60 minutes each time. Continue treatment until recovery.

39. Constipation:

Irradiate **Tianshu acupoint(Pic 51)**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

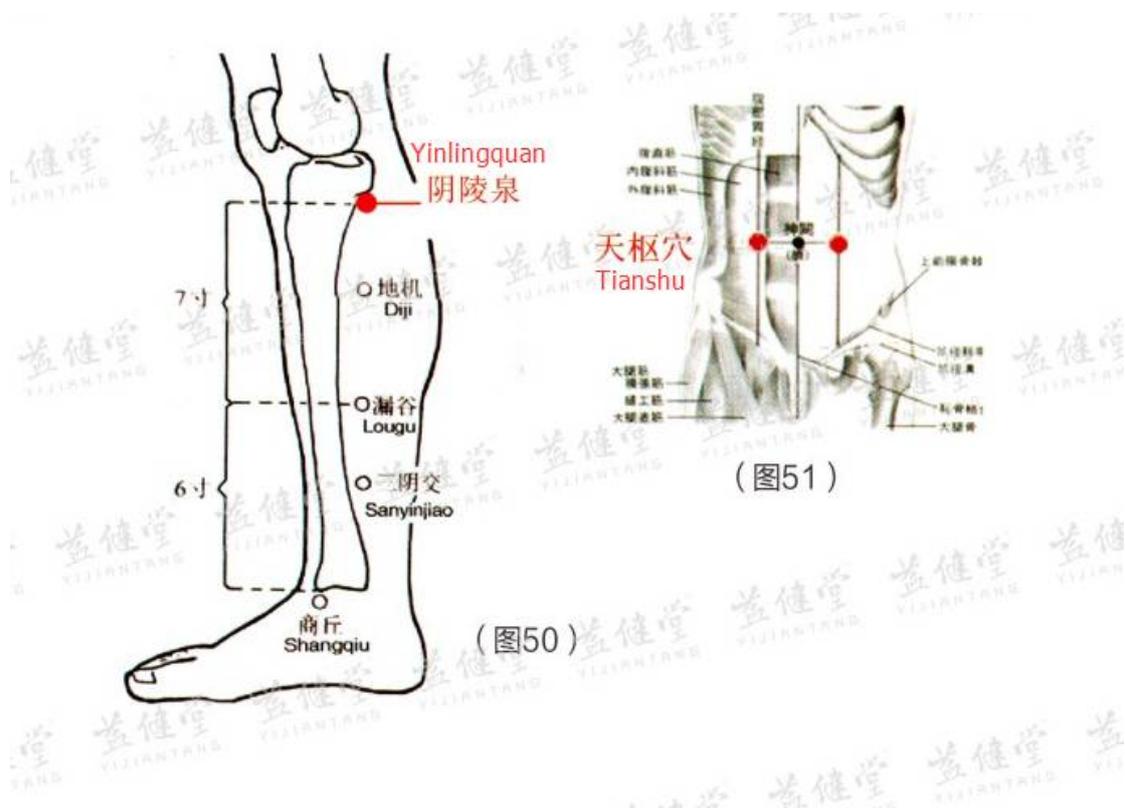

40. Insomnia, neurasthenia

Irradiate **Anmian acupoint (Pic 52)**, **Neiguan acupoint (Pic 53)**. 1 or 2 times a day, 30 minutes each time. Continue treatment until recovery.

41. Dysmenorrhea:

Irradiate underbelly. 1 or 2 times a day, 30 minutes each tim. Continue treatment until recovery.



(图52)

(图53)

42. Pelvic inflammation, Pelvic inflammatory mass, Accessory phlogistic, Senile vaginitis:

Irradiate the pain points of **hypogastrium and Huiyin(Pic 54)** Acupuncture. 1 or 2 times a day, 30 minutes each time, continued treatment until recovery.

Huiyin Acupuncture: located in the middle sag of the anus and genitals.

43. Acne:

Irradiate the acne area. 1 or 2 times a day, 30 minutes each time, continued use until recovery.


44. Eczema, Herpes zoster, Herpes zoster neuralgia, Allergic dermatitis

Irradiate the skin lesion area.

1 or 2 times a day, 30 minutes each time, continued use until recovery.





(图54)

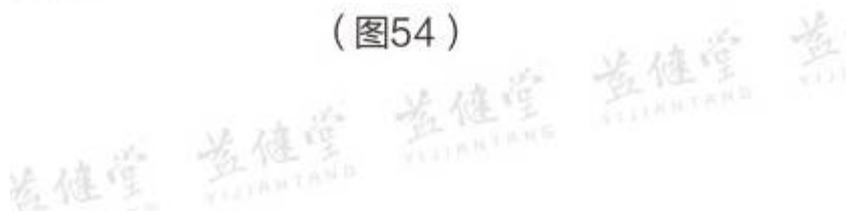

II. Millimeter-wave cancer treatment


1. Cancer of head and neck:

Irradiate **Hegu, Yingxiang, Xiache and Ashi**(No fixed position, around the affected area, pressing pain point) accupoint. One time a day, 30 minutes each time, 5 times a week, 30 times a course.

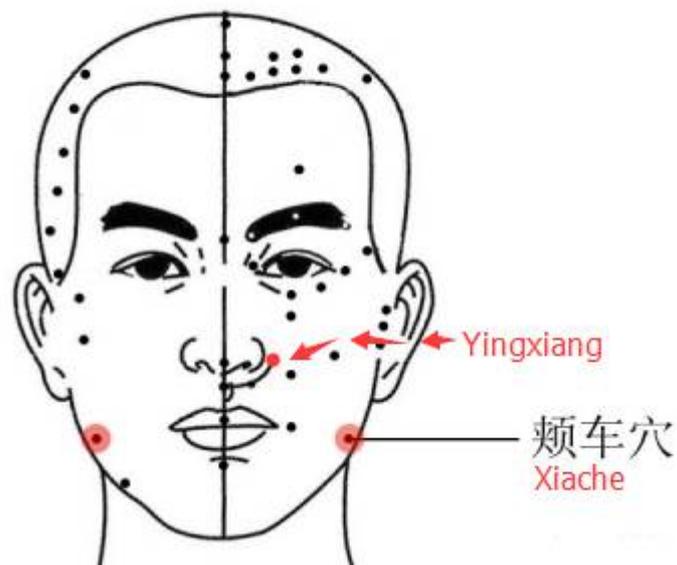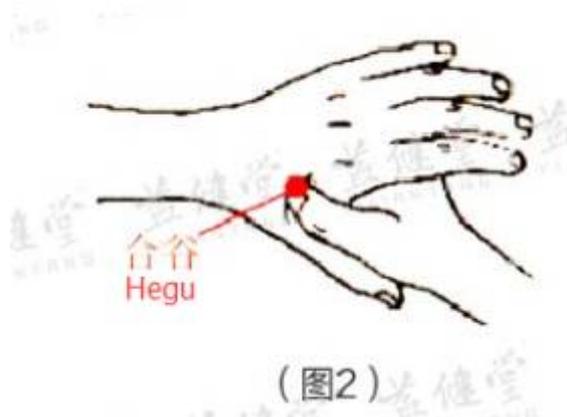

2. Breast cancer:

Irradiate **Rugen**, **Danzhong**, **Neiguan** accupoint. One time a day, 30 minutes each time, 5 times a week, 30 times a course.

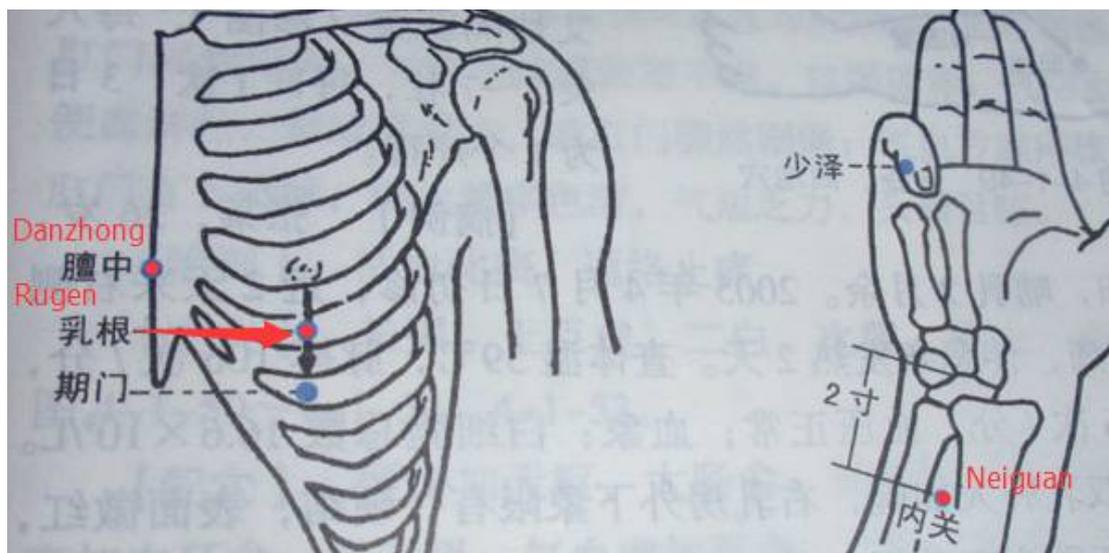


3. Lung cancer:

Irradiate **Tiantu**, **Feishu**, **Danzhong** accupoint. One time a day, 30 minutes each time, 5 times a week, 30 times a course.




4. Rectal cancer:


Irradiate **Qihai, Guanyuan and Zusanli** accupoint. One time a day, 30 minutes each time, 5 times a week, 30 times a course.

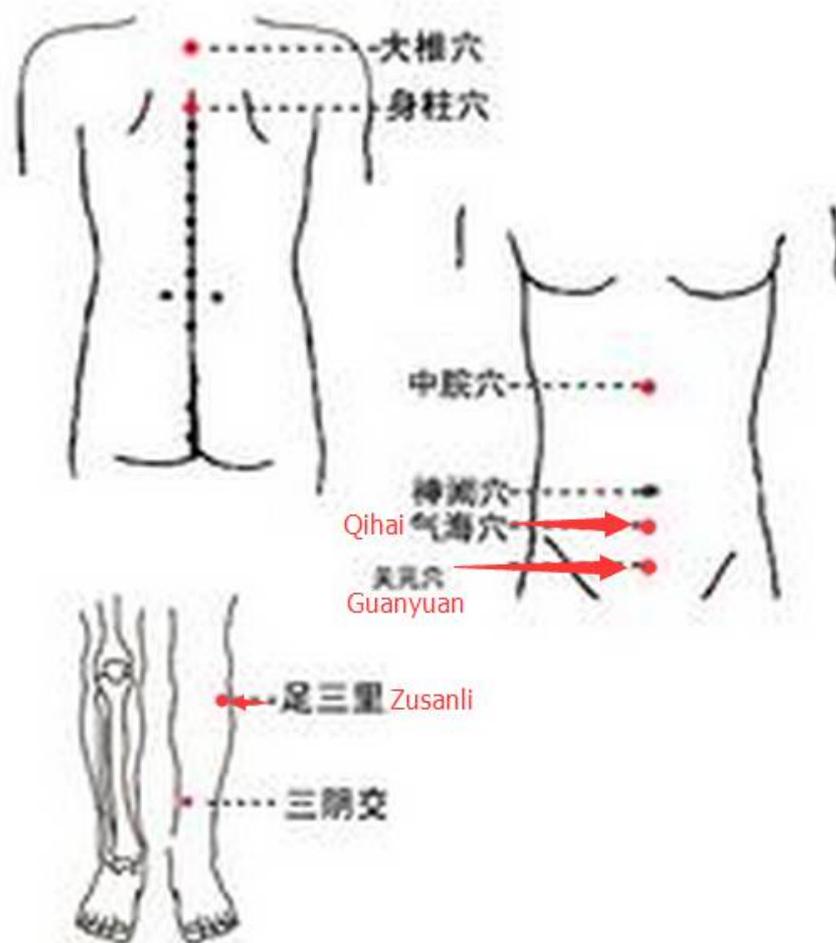

5. Bladder cancer:

Irradiate **Zhongji, Panguangshu, Yinlingquan** accupoint. One time a day, 30 minutes each time, 5 times a week, 30 times a course.

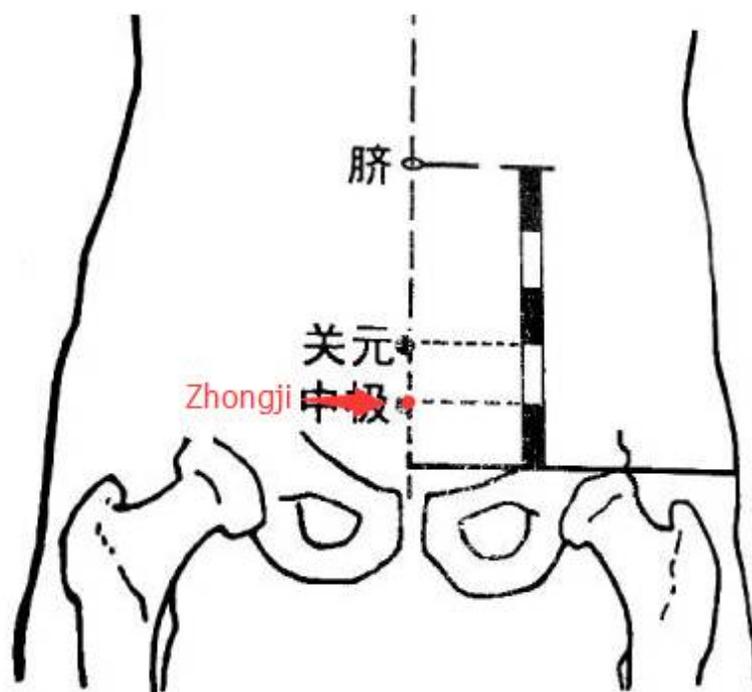

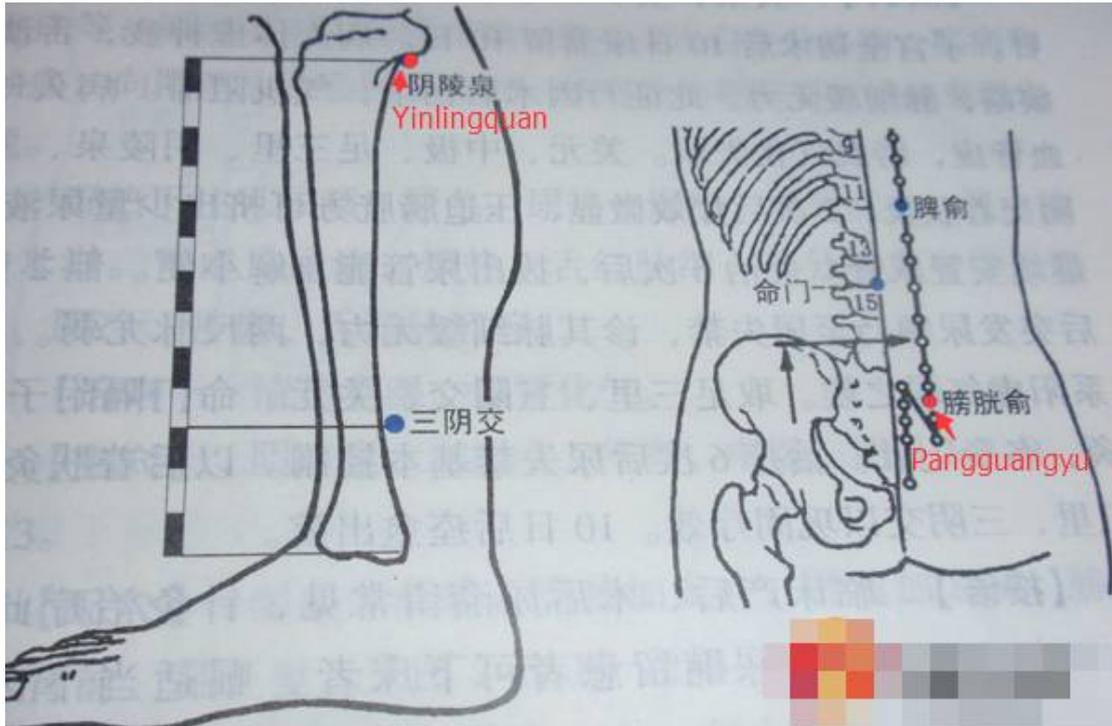

International School of Healing Arts

In Möslim, last evening prayer - In Hindu, Goddess of protection - In Arabic, a long life

Reserving all rights, defences and typemistakes.

Yours sincerely, ISHA Global Sales Team,

Email: sales@nonlinearsystem.eu

- **Headquarters Office The Netherlands:**
ISHA Quantum Metaphysics of NLS
The Netherlands, Europe
- **Email at** support@nonlinearsystem.eu or sales@nonlinearsystem.eu
- **TICKET SUPPORT** <https://nonlinearsystem.eu/support/>

Disclaimer: The information is provided for educational, experimental, investigational, and research purposes only.

Disclaimer: Information contained on this email should not be construed as medical advice.

By ordering you accept our terms and conditions <http://nonlinearsystem.eu/terms-and-conditions> and our medical disclaimer <https://nonlinearsystem.eu/medical-disclaimer>
----- end -----